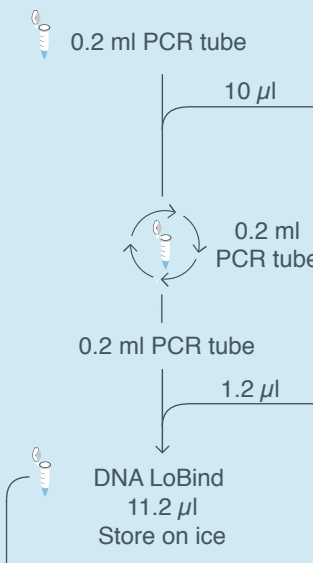


Rapid sequencing of genomic DNA for the MinION™ device using SQK-RAD002 (1/2)

Flow Cell Number
DNA Samples


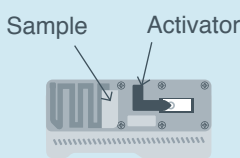

Before start checklist

- | | | |
|---|---|--|
| <input type="checkbox"/> Rapid Sequencing Kit (SQK-RAD002 with EXP-LLB001) | <input type="checkbox"/> 0.2 ml thin-walled PCR tubes | <input type="checkbox"/> Timer |
| <input type="checkbox"/> Pipettes and tips P1000, P200, P20, P10 and P2 1.5ml | <input type="checkbox"/> MinION SpotON Flow cell (FLO-MIN106) | <input type="checkbox"/> Microfuge |
| <input type="checkbox"/> Eppendorf DNA LoBind tubes | <input type="checkbox"/> Nuclease-free water (NFW) | <input type="checkbox"/> Thermal cycler at 30 °C and 75 °C |
| | <input type="checkbox"/> NEB Blunt/TA Ligase Master Mix (MO367) | |

MASSFLOW	INSTRUCTIONS	NOTES / OBSERVATIONS
	<input type="checkbox"/> Take 200 ng high molecular weight DNA in 7.5 µl <input type="checkbox"/> Add 2.5 µl FRM Mix gently by inversion + spin down Incubate for 1 min at 30 °C then 1 min at 75 °C Spin down briefly	FRM
	<input type="checkbox"/> Add 1 µl RAD <input type="checkbox"/> Add 0.2 µl Blunt/TA Ligase Master Mix Incubate for 5 mins at RT The library preparation is complete and ready for loading onto the MinION	RAD

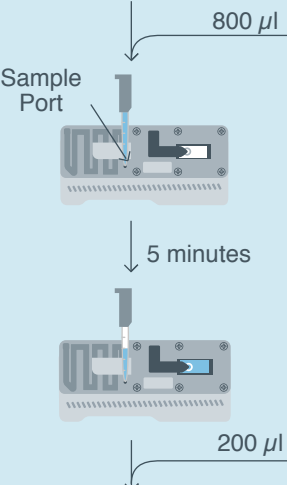
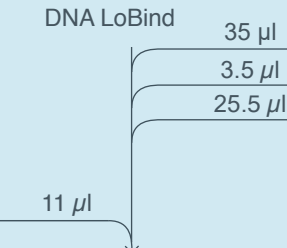
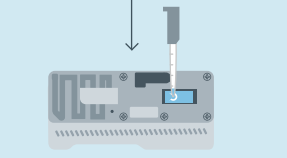
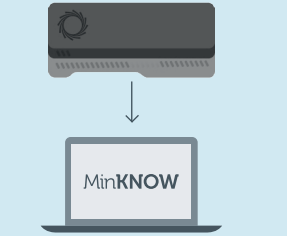
Before start checklist

- | | | |
|--|--|--|
| <input type="checkbox"/> MinION™ connected to computer with SpotON Flow Cell | <input type="checkbox"/> Computer set up to run MinKNOW | <input type="checkbox"/> PSM, RBF and LLB on ice |
| <input type="checkbox"/> Run platform QC in parallel to library prep | <input type="checkbox"/> Desktop Agent set up | <input type="checkbox"/> NFW at RT |
| | <input type="checkbox"/> Run Name set | <input type="checkbox"/> Platform QC completed |
| | <input type="checkbox"/> Pre-sequencing Mix (PSM) at > 4 ng/µl | |

Priming and loading the library   	Prepare the MinION for sequencing protocol This step can be run in parallel with the preparation of the library from genomic DNA to Pre-sequencing Mix <input type="checkbox"/> Assemble the MinION and MinION Flow Cell <input type="checkbox"/> Setup MinKNOW to run the Platform QC – name the run and start the protocol script – NC_Platform_QC.py <input type="checkbox"/> Allow the script to run to completion and the number of active pores are reported	
	Prime the Flow Cell ready for the library to be loaded when library preparation is complete Prepare priming buffer <input type="checkbox"/> 480 µl RBF <input type="checkbox"/> 520 µl Nuclease-free water	RBF

Rapid sequencing of genomic DNA for the MinION™ device using SQK-RAD002 (2/2)

Flow Cell Number
DNA Samples

MASSFLOW	INSTRUCTIONS	NOTES / OBSERVATIONS
 <p>800 µl</p> <p>Sample Port</p> <p>5 minutes</p> <p>200 µl</p>	<p>Prime the Flow Cell</p> <p>Open the sample port. Draw back a few µls of buffer to make sure there is continuous buffer flow from the sample port across the sensor array.</p> <ul style="list-style-type: none"> <input type="checkbox"/> Load 800 µl of the priming buffer. Wait 5 minutes <input type="checkbox"/> Gently lift the activator to make the SpotON port accessible <input type="checkbox"/> Load 200 µl of the priming buffer as before 	
 <p>DNA LoBind 35 µl</p> <p>3.5 µl</p> <p>26.5 µl</p> <p>11 µl</p>	<p>Prepare the library for loading</p> <ul style="list-style-type: none"> <input type="checkbox"/> 25.5 µl RBF kept on ice <input type="checkbox"/> 12 µl NFW kept at RT <input type="checkbox"/> 26.5 µl LLB kept on ice <input type="checkbox"/> 11 µl Adapted and tethered library <p>Mix by inversion and spin down</p>	<div> <div>RBF</div> <div>LLB</div> </div>
 <p>DNA LoBind 75 µl</p>	<p>Loading the prepared library</p> <ul style="list-style-type: none"> <input type="checkbox"/> Add 75 µl of sample to the flow cell via the SpotON port in a dropwise fashion. Ensure each drop flows into the port before adding the next. <input type="checkbox"/> Gently replace the activator, making sure the bung enters the SpotON port <input type="checkbox"/> Close the sample port cover and replace the MinION lid. 	
 <p>MinKNOW</p>	<p>Starting the sequencing script in MinKNOW</p> <ul style="list-style-type: none"> <input type="checkbox"/> Return to MinKNOW, name the run, select the NC_48Hr_Sequencing_FLO_MIN106_SQKLSK108_plus_Basecaller.py for live basecalling using the start in the MinKNOW dialogue box <input type="checkbox"/> MinKNOW will report the number of pores available for sequencing before data collection begins. These may differ from those reported in the Platform QC. <input type="checkbox"/> Allow the protocol to proceed until MinKNOW reports Finished Successfully. Use the Stop in the Control Panel to finish the protocol. <input type="checkbox"/> Close down MinKNOW and disconnect the MinION <p>If using Albacore for local basecalling please refer to the instructions in Albacore basecalling software</p>	

Before start checklist

- ☐ Store washed flow cell at 4 °C or complete the returns form in the library prep
- ☐ Store MinION at RT
- ☐ Return reagents to the freezer
- ☐ Navigate to www.metrichor.com to review the full sequencing report